A E N O

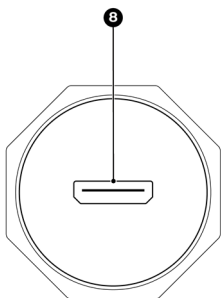
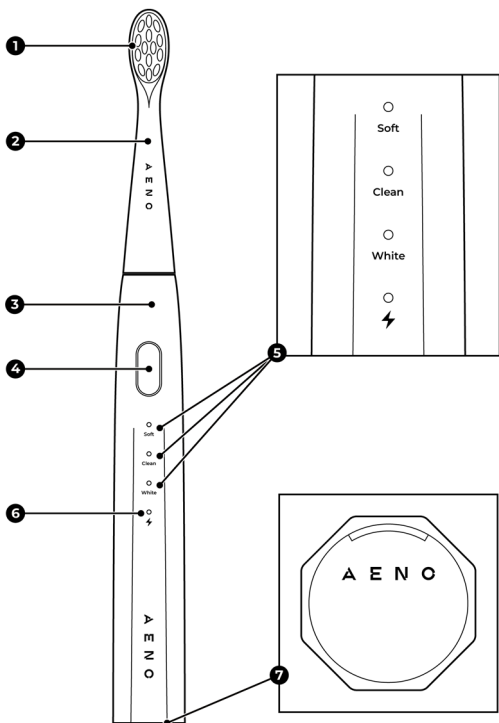
Quick Start Guide

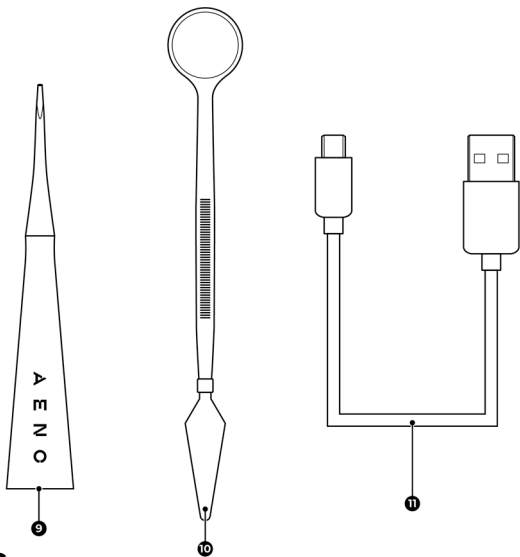
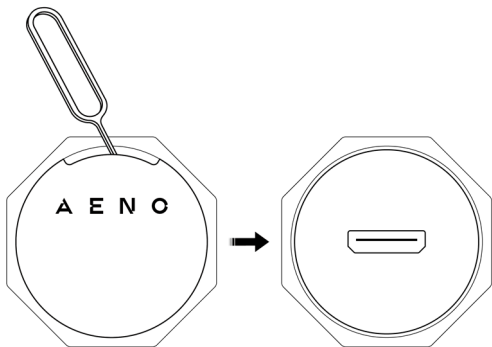
Toothbrush DB7 / DB8

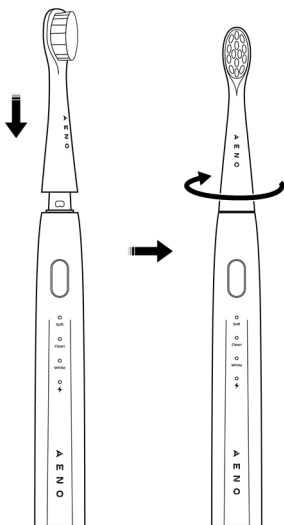
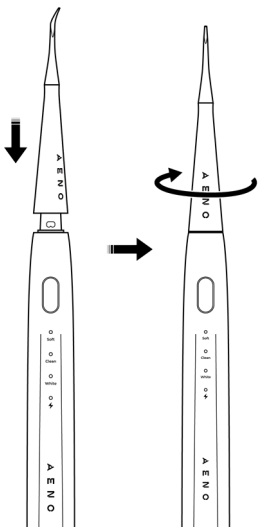
(Models: ADB0007,
ADB0008)

12.01.2024 | v 2.0.1

| | |
|----------|----|
| ENG..... | 6 |
| ARA..... | 8 |
| BOS..... | 10 |
| BUL..... | 12 |
| CES..... | 14 |
| DEU..... | 16 |
| ELL..... | 18 |
| EST..... | 20 |
| FRA..... | 22 |
| HRV..... | 24 |
| HUN..... | 26 |
| HYE..... | 28 |
| ITA..... | 30 |
| KAT..... | 32 |
| KAZ..... | 34 |
| LAV..... | 36 |
| LIT..... | 38 |
| NLD..... | 40 |
| POL..... | 42 |
| POR..... | 44 |
| RON..... | 46 |
| RUS..... | 48 |
| SLK..... | 50 |
| SLV..... | 52 |
| SPA..... | 54 |
| SRP..... | 56 |
| UKR..... | 58 |
| UZB..... | 60 |

B**C**

D**E**

ENG AENO™ Toothbrush is intended for daily oral hygiene.

Technical Specifications

Power supply (input): 5.0 V (DC), 0.3 A. Power: 1.5 W (max.). Battery: non-removable, lithium-ion; 800 mAh, 3.7 V. Battery full charge time: up to 4 h. Operating time on a single charge: up to 100 days. Ingress protection rating: IPX7. Noise level (at a distance of 0.3 m): ≤ 60 dB. Number of operating modes: 3. Vibration intensity: up to 30,000 vibrations per minute. Colour: white. Dimensions of the toothbrush with the bristle brush head (diameter × height): 20.8×208 mm. USB (A) / micro-USB cable length: 0.3 m. Operating conditions: temp. 0...+35 °C, RH < 60 %. Storage conditions: temp. 0...+45 °C, RH < 60 %.

Scope of Delivery (see fig. A, B)

ADB0007: Toothbrush (handle + replacement head with bristles (1 pc.)), USB (A) / micro-USB cable, sticker sheets (2 pcs.), Quick Start Guide.

ADB0008: Toothbrush (handle + replacement heads with bristles (3 pcs.)), plaque remover, dental mirror, USB (A) / micro-USB cable, Quick Start Guide.

Device Elements and Accessories (see fig. A, B)

1 – bristles of the replacement brush head, 2 – replacement head with bristles, 3 – handle, 4 – power button, 5 – operating mode indicators, 6 – battery charge indicator, 7 – protective cover, 8 – micro-USB port, 9 – plaque remover, 10 – dental mirror, 11 – USB (A) / micro-USB cable.

Restrictions and Warnings

This device is not intended for use by children under the age of 6 or by persons with reduced physical, mental or intellectual abilities unless they have sufficient experience and knowledge of how to use the device and are supervised by a person responsible for their safety. Do not allow children to play with the device and its accessories. Do not drop or throw the device. Do not turn on or use a damaged device. If you have a pacemaker, oral diseases, or if you have just undergone oral surgery, consult your doctor before using the device. If your gums begin to bleed after using the device and the bleeding continues for several days, contact your dentist. Each family member should use an individual brush head. Replace the brush head at least once every 3 months. Use only original replacement heads. The tartar remover is intended for use by adults only. If you are not using the device for a long time, charge it every 2–3 months to prevent deep discharge of the battery. Do not wash the device or its accessories in a dishwasher. Do not use chemical and aggressive detergents, abrasive pastes, products containing acids and solvents or metal sponges for cleaning.

Statuses of the Battery Charge Indicator

| Colour | Status | Meaning |
|--------|--------------|---|
| Green | Steadily lit | The device is turned on. Battery level is medium/high |
| | | The device is connected to a power source. Battery is fully charged |
| Yellow | Steadily lit | The device is connected to a power source. Charging is in progress |
| Red | Steadily lit | The device is turned on. Battery level is low |
| | Blinking | The battery is fully discharged |
| - | Off | The device is not connected to a power source or is turned off |

Statuses of the Operating Mode Indicators

You can tell the current operating mode of the toothbrush by the steadily lit white indicator located above the name of the mode.

Power Button Functions

| Action | Meaning |
|--|------------------------|
| Press once (device is turned off) | Turn the device on |
| Press once (the device is turned on, and it has been less than 5 seconds since it was turned on) | Switch operating modes |
| Press once (the device is turned on, and it has been more than 5 seconds since it was turned on) | Turn the device off |

Preparing for Operation

Carefully open the packaging and remove the device and its accessories. Clean them (see "Cleaning and Maintenance"). Remove the protective cover from the brush handle (see fig. C).

WARNING! The protective cover ensures that the micro-USB port is waterproof. If the protective cover is removed, do not allow liquids to enter the port of the device.

Connect the device to a power source using a USB (A) / micro-USB cable. The battery charge indicator will light up yellow. When the indicator turns green, the battery is fully charged. Disconnect the device from the power source, and then disconnect the cable from the micro-USB port. Insert the protective cover into the device handle.

Choosing the Operation Modes of the Device

Turn on the device by pressing the power button once. Subsequent presses of the button will switch the operating modes. The last used mode is saved in the device settings the next time the device is turned on.

Operating Modes

- "Soft". Delicate, gentle brushing of the teeth. Recommended for brushing sensitive teeth, as well as for brushing children's teeth.
- "Clean". Daily brushing. Recommended for standard brushing.
- "White". Teeth whitening. Recommended for people who smoke and drink tea or coffee often.

Operation

With the head with bristles. Install the head as shown in fig. **D**. Wet the bristles of the brush head and apply the desired amount of a dental care product to the bristles. Point the bristles toward the surface of the teeth at an angle of about 45°. Press the power button to turn on the device. Brush your teeth without putting too much pressure on the enamel. Every 30 seconds, the device will notify you to change the cleaning zone by briefly stopping the vibration. Change the cleaning zone after each notification. After 2 minutes of operation, the device will turn off automatically.

With plaque remover. Install the attachment as shown in fig. **E**. Point the metal end of the attachment toward the tooth surface at an angle of about 45°. Press the power button to turn on the device. Clean problem areas with gentle vertical strokes.

WARNING! Do not apply excessive force when using the plaque remover. Careless use of the attachment can cause oral trauma as well as damage tooth enamel.

Every 30 seconds, the device will notify you to change the cleaning zone by briefly stopping the vibration. Change the cleaning zone after each notification. After 2 minutes of operation, the device will turn off automatically.

Note. Using a plaque remover cannot substitute professional oral hygiene at the dentist.

Cleaning and Maintenance

WARNING! Do not clean the device until it has been turned off and disconnected from the power source.

Rinse the device and accessories under running water. Then wipe dry.

Troubleshooting

1. The toothbrush does not turn on. Possible cause: the battery is discharged. Solution: charge the device (charging time – at least 4 hours).
2. The device is connected to a power source but is not charging (the battery charge indicator is not lit). Possible cause: the device has not been used for a long time. Solution: connect the toothbrush to a power source and wait about 15 minutes. If the battery charge indicator does not light up, contact an authorized service centre.

If you have any questions or difficulties in using the AENO™ device, please e-mail support at support@aeno.com or chat with us online at aeno.com/service-and-warranty. The specialists will help you solve the problem, and you will not have to waste time and effort visiting the store.

ASBISc reserves the right to modify the device and to make changes and additions to this document without prior notice to users.

The warranty period and service life is 2 years from the date of retail sale of the device.

Manufacturer: ASBISc Enterprises PLC, Iapetou, 1, Agios Athanasios, 4101, Limassol, Cyprus. Made in China. All trademarks and brands mentioned herein are the property of their respective owners. The date of manufacture is indicated on the package.

Up-to-date information and a detailed description of the device, as well as connection instructions, certificates, information on companies that accept quality claims and warranties are available at aeno.com/documents.

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ENG Recycling information These symbols indicate that you must follow the Waste Electrical and Electronic Equipment (WEEE) and Waste Battery and Battery Regulations when disposing of the device, its batteries and accumulators, and its electrical and electronic accessories. According to the rules, this equipment must be disposed of separately at the end of its service life. Do not dispose of the device, its batteries and accumulators, or its electrical and electronic accessories together with unsorted municipal waste, as this will harm the environment. To dispose of this equipment, it must be returned to the point of sale or turned in to a local recycling facility. You should contact your local household waste disposal service for details.

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