aquapick[®]

SONIC ELECTRIC TOOTHBRUSH

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Advanced Sonic Electric Toothbrush
Periuța Sonică Electrică Aquapick

IMPORTANT SAFEGUARDS

The Sonic electric toothbrush is safe and effective when used according to directions.

READ ALL INSTRUCTIONS BEFORE USING:

1. DANGER-WATER

- 1 Do not place or store the charger base where it can fall or be pulled into a bathtub, sink or toilet.
- 2 Do not place or drop the charger into water or other liquid.
- 3 Do not use while bathing.
- 4 Unplug charger base before cleaning and assure that it is dry before plugging it into an usb port.

2. DANGER-WARNINGS

- 1 Use this product only for its intended use as described in this manual.
- 2 This product is designed to clean your teeth, gums and tongue only.
Do not use any attachments or brush heads other than those recommended by the manufacturer.
- 3 If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean your brush head and handle to avoid cracking the plastic.
- 4 Never operate this product if it has a damaged cord, plug or charger base. Do not use this product if it is not working properly, if it has been dropped or damaged or dropped into water. Immediately discontinue using if the product appears damaged in any way.
- 5 Keep the toothbrush, charger base and power cord away from heated surfaces.
- 6 Do not attempt to remove the battery.

3. MEDICAL-WARNINGS

- 1 If you have had oral surgery or gum surgery within the past two months, consult your dental provider before using the product.
- 2 Consult your dental provider if excessive bleeding occurs while using this product or minor bleeding continues to occur after one week of use.
- 3 Consult your physician or implant device manufacturer prior to using this product if you have questions concerning the combined use of this product with a pacemaker or other implanted medical devices.
- 4 Consult your physician prior to using the toothbrush if you have any medical concerns.
- 5 Never use while sleeping or drowsy.

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

CONFIGURATION

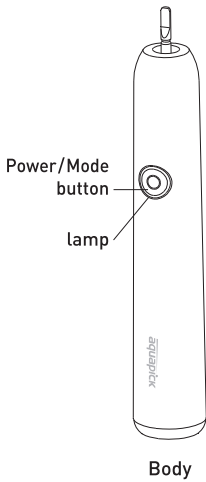
* Purchasable Separately
1 SET (4pcs/set)

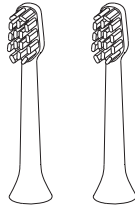Brush head (2pcs)

User manual

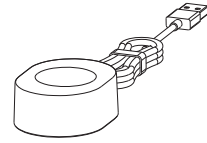Charging stand

GETTING STARTED

Charging the toothbrush

When the product is being charged, the charging indicator flashes slowly in white and glows blue when fully charged.

Notes:

- 1 It takes 12 hours to charge the battery fully. (Fig.1)
- 2 The toothbrush must be charged for 12 hours, if:
 - a. this is the first time to use the toothbrush;
 - b. it has not been used for more than three months;
- 3 We recommend that you keep the toothbrush in the charger when not in use.
- 4 The operating duration will become shorter as the battery gets older. If the toothbrush operates only for short time after being fully charged, this means that the battery is nearing the end of its life cycle.

Contact service center for support.




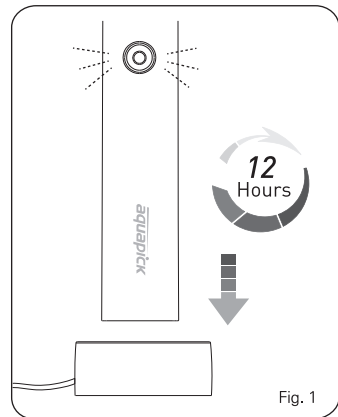
Fig. 1

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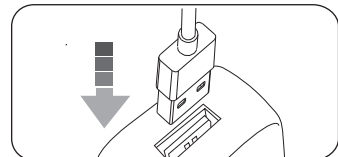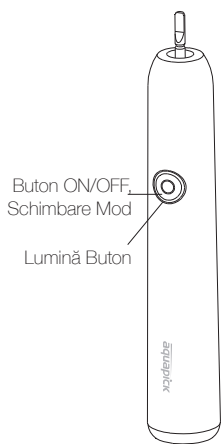
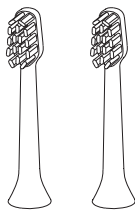
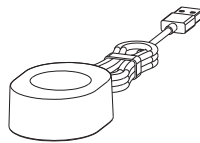
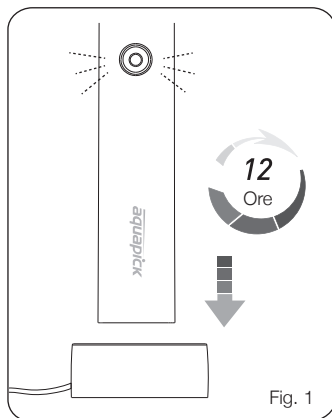
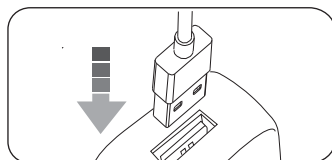
HOW TO USE


- 1 Slide the brush head into the handle.
Push brush head down until it locks into place on handle.
Note: There is small gap between the brush head and the handle.(Fig. 2)

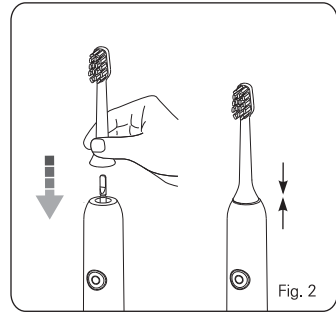- 2 Wet the bristles and apply a small amount of toothpaste.
- 3 Place the bristles against the teeth at a slight angle towards the gum line.(Fig. 3)

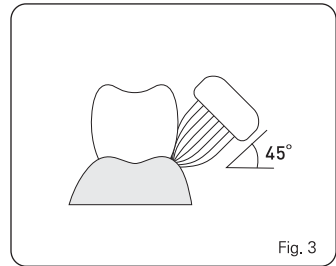
- 4 The brushing modes are as follows:
 - a. Clean Mode: Standard mode for superior teeth cleaning
 - b. Sensitive Mode: Gently cleans your teeth or any other sensitive areas;
- 5 If the toothbrush runs under "Clean" mode, under which the light ring illuminates in white, and the running time is less than 5 seconds, press the Switch/Function button once, the toothbrush will switch to "Sensitive" mode and the light ring will illuminate in blue. If the running time exceeds 5 seconds, press the Switch/Function button again, the toothbrush will turn off and the current working mode will be saved, which means that the toothbrush will work under "Clean" mode after being activated next time.

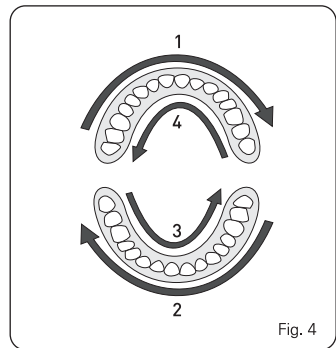- 6 If the toothbrush runs under "Sensitive" mode, under which the light ring illuminates in blue, and the running time is less than 5 seconds, press the Switch/Function button once, the toothbrush will switch to "Clean" mode and the light ring will illuminate in white. If the running time exceeds 5 seconds, press the Switch/Function button once, the toothbrush will turn off and the current working mode will be saved, which means that the toothbrush will work under "Sensitive" mode after being activated next time.

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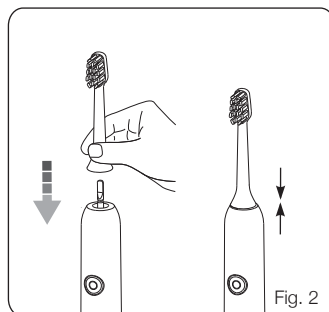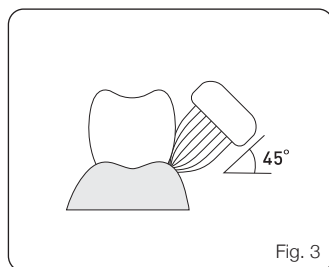
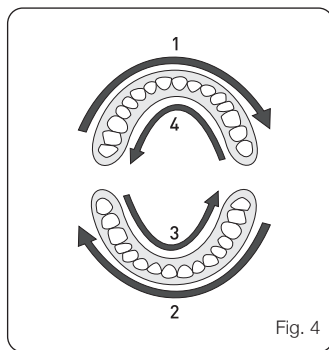
- 7 Apply light pressure and gently move the brush head slowly across the teeth in a small back and forth motion.
- 8 In order to maximize cleaning effectiveness, divide your mouth in 4 sections. Start with Section 1 and brush each section for 30 seconds. (Fig.3) The brushing action will pause between each 30 second interval, indicating when to move on to the next section. The toothbrush will automatically stop after 2 minutes of use.
- 9 Pull brush head out from the handle. (Fig 5)
- 10 Rinse both the brush head and handle in running water and dry them. (Fig 6)
- 11 Put the toothbrush back into the charger.
- ※ It is better to sterilize brush heads by putting them in Aquapick toothbrush sterilizer, (separate sale)

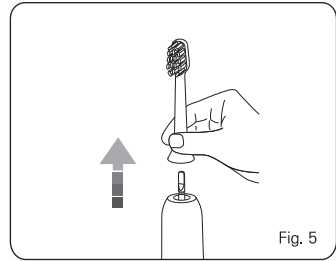
Fig. 5



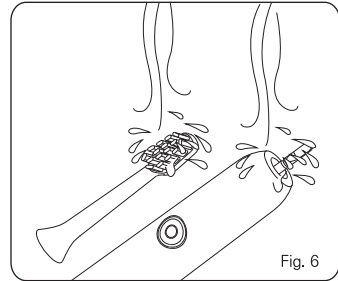
Fig. 6

MAINTENANCE

Toothbrush

- 1 Remove the brush head.
- 2 Wipe all surface areas with a damp cloth and put the toothbrush back into the charger.
 - ※ Do not submerge the toothbrush into water.
- 3 Rinse the brush head and bristles after each use.
- 4 Do not clean in the dishwasher.

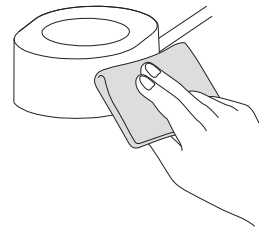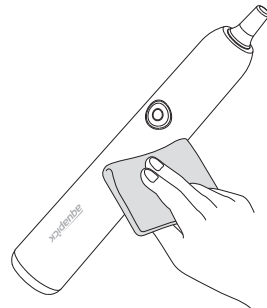Charger

Using a soft cloth, wipe off any dirt on the charger and the storage.

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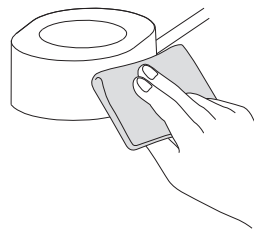
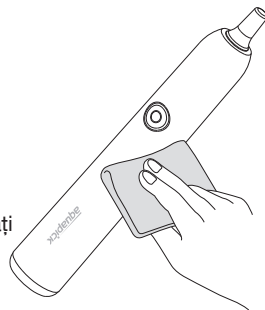
TROUBLE SHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
The toothbrush doesn't work.	You have just purchased the product or it has not been used for more than 3 months.	Charge the main unit for at least 12 hours.
Weak or diminished vibration.	The toothbrush is pressed hard against your teeth.	Place the brush head gently against your teeth.
The toothbrush operates only for a short period of time even immediately after being charged.	The battery has reached the end of its life cycle.	Contact service center for support.
	Charging time was insufficient.	Charge the main unit for at least 12 hours.

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aquapick[®]
SONIC ELECTRIC TOOTHBRUSH